## ORGENTEC Diagnostika GmbH

Carl-Zeiss-Straße 49-51

55129 Mainz - Germany

Phone: +49 (0) 61 31 / 92 58-0

Fax: +49 (0) 61 31 / 92 58-58

Internet: www.orgentec.com

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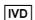
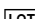
2015-02

## ORG 243M Anti-Annexin V IgM

### KRÁTKÝ POPIS

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## REFERENCE

- Banzato A, Pozzi N, Frasson R, De F, V, Ruffatti A, Bison E et al. Antibodies to Domain I of beta(2)Glycoprotein I are in close relation to patients risk categories in Antiphospholipid Syndrome (APS). *Thromb Res* 2011; 128 (6):583-6.
- Bertolaccini ML, Amengual O, Atsumi T, Binder WL, de LB, Forastiero R et al. 'Non-criteria' aPL tests: report of a task force and preconference workshop at the 13th International Congress on Antiphospholipid Antibodies, Galveston, TX, USA, April 2010. *Lupus* 2011; 20(2):191-205.
- de Laat B, de Groot PG. Autoantibodies directed against domain I of beta2-glycoprotein I. *Curr Rheumatol Rep* 2011; 13(1):70-6.
- de Laat B, Mertens K, de Groot PG. Mechanisms of disease: antiphospholipid antibodies-from clinical association to pathologic mechanism. *Nat Clin Pract Rheumatol* 2008; 4(4):192-9.
- de Laat B, Pengo V, Pabinger I, Musial J, Voskuyl AE, Bultink IE et al. The association between circulating antibodies against domain I of beta2-glycoprotein I and thrombosis: an international multicenter study. *J Thromb Haemost* 2009; 7(11):1767-73.
- Espinosa G, Cervera R. Antiphospholipid syndrome. *Arthritis Res Ther* 2008; 10(6):230.
- Favaloro EJ, Wong RC. Laboratory testing for the antiphospholipid syndrome: making sense of antiphospholipid antibody assays. *Clin Chem Lab Med* 2011; 49(3):447-61.
- Fischer MJ, Rauch J, Levine JS. The antiphospholipid syndrome. *Arthritis Rheum* 2007; 27(1):35-46.
- Giannakopoulos B, Passam F, Ioannou Y, Krilis SA. How we diagnose the antiphospholipid syndrome. *Blood* 2009; 113(5):985-94.
- Greaves M, Cohen H, Machin SJ, Mackie I. Guidelines on the investigation and management of the antiphospholipid syndrome. *Br J Haematol* 2000; 109(4):704-15.
- Hughes GR. Hughes syndrome: antiphospholipid syndrome. *J R Coll Physicians Lond* 1998; 32(3):260-4.
- Hughes GR. Hughes Syndrome (the antiphospholipid syndrome): ten clinical lessons. *Autoimmun Rev* 2008; 7 (3):262-6.
- Hughes GR. Antiphospholipid syndrome, migraine and stroke. *Lupus* 2010; 19(5):555-6.
- Hughes GR, Harris NN, Gharavi AE. The anticardiolipin syndrome. *J Rheumatol* 1986; 13(3):486-9.
- Koike T, Bohgaki M, Amengual O, Atsumi T. Antiphospholipid antibodies: lessons from the bench. *J Autoimmun* 2007; 28(2-3):129-33.
- Lakos G, Favaloro EJ, Harris EN, Meroni PL, Tincani A, Wong RC et al. International consensus guidelines on anticardiolipin and anti-beta2-glycoprotein I testing: report from the 13th International Congress on Antiphospholipid Antibodies. *Arthritis Rheum* 2012; 64(1):1-10.
- Mackworth-Young C. Primary antiphospholipid syndrome: a distinct entity? *Autoimmun Rev* 2006; 5(1):70-5.
- Miyakis S, Lockshin MD, Atsumi T, Branch DW, Brey RL, Cervera R et al. International consensus statement on an update of the classification criteria for definite antiphospholipid syndrome (APS). *J Thromb Haemost* 2006; 4(2):295-306.
- Molina JF, Gutierrez-Urena S, Molina J, Uribe O, Richards S, De CC et al. Variability of anticardiolipin antibody isotype distribution in 3 geographic populations of patients with systemic lupus erythematosus. *J Rheumatol* 1997; 24(2):291-6.
- Oku K, Atsumi T, Amengual O, Koike T. Antiprothrombin antibody testing: detection and clinical utility. *Semin Thromb Hemost* 2008; 34(4):335-9.
- Pengo V, Biasiolo A, Bison E, Chantarangkul V, Tripodi A. Antiphospholipid antibody ELISAs: survey on the performance of clinical laboratories assessed by using lyophilized affinity-purified IgG with anticardiolipin and anti-beta2-Glycoprotein I activity. *Thromb Res* 2007; 120(1):127-33.
- Pierangeli SS, de Groot PG, Dlott J, Favaloro E, Harris EN, Lakos G et al. 'Criteria' aPL tests: report of a task force and preconference workshop at the 13th International Congress on Antiphospholipid Antibodies, Galveston, Texas, April 2010. *Lupus* 2011; 20(2):182-90.

23. Pierangeli SS, Favaloro EJ, Lakos G, Meroni PL, Tincani A, Wong RC et al. Standards and reference materials for the anticardiolipin and anti-beta-2-glycoprotein I assays: a report of recommendations from the APL Task Force at the 13th International Congress on Antiphospholipid Antibodies. Clin Chim Acta 2012; 413(1-2):358-60.
24. Sinico RA, Bollini B, Sabadini E, Di Toma L, Radice A. The use of laboratory tests in diagnosis and monitoring of systemic lupus erythematosus. J Nephrol JID - 9012268 2002; 15 Suppl 6:S20-S27.
25. Tincani A, Andreoli L, Casu C, Cattaneo R, Meroni P. Antiphospholipid antibody profile: implications for the evaluation and management of patients. Lupus 2010; 19(4):432-5.
26. Tincani A, Morozzi G, Afeltra A, Alessandri C, Allegri F, Bistoni O et al. Antiprothrombin antibodies: a comparative analysis of homemade and commercial methods. A collaborative study by the Forum Interdisciplinare per la Ricerca nelle Malattie Autoimmuni (FIRMA). Clin Exp Rheumatol 2007; 25(2):268-74.
27. Wilson WA, Gharavi AE, Koike T, Lockshin MD, Branch DW, Piette JC et al. International consensus statement on preliminary classification criteria for definite antiphospholipid syndrome: report of an international workshop. Arthritis Rheum 1999; 42(7):1309-11.
28. Wong RC, Favaloro EJ, Adelstein S, Baumgart K, Bird R, Brighton TA et al. Consensus guidelines on anti-beta 2 glycoprotein I testing and reporting. Pathology 2008; 40(1):58-63.
29. Wong RC, Gillis D, Adelstein S, Baumgart K, Favaloro EJ, Hendle MJ et al. Consensus guidelines on anticardiolipin antibody testing and reporting. Pathology 2004; 36(1):63-8.