## ORGENTEC Diagnostika GmbH

Carl-Zeiss-Straße 49-51

55129 Mainz - Germany

Phone: +49 (0) 61 31 / 92 58-0

Fax: +49 (0) 61 31 / 92 58-58

Internet: www.orgentec.com

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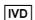
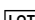
2014-06

## ORG 232 Anti-RNP-70

### KRÁTKÝ POPIS

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
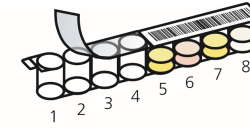
### REFERENCE

- Alba P, Bento L, Cuadrado MJ, Karim Y, Tungekar MF, Abbs I et al. Anti-dsDNA, anti-Sm antibodies, and the lupus anticoagulant: significant factors associated with lupus nephritis. *Ann Rheum Dis* 2003; 62(6):556-560.
- Antico A, Platzgummer S, Bassetti D, Bizzaro N, Tozzoli R, Villalta D. Diagnosing systemic lupus erythematosus: new-generation immunoassays for measurement of anti-dsDNA antibodies are an effective alternative to the Farr technique and the Crithidia luciliae immunofluorescence test. *Lupus* 2010; 19(8):906-912.
- Brouwer R, Hengstman GJ, Vree EW, Ehrfeld H, Bozic B, Ghirardello A et al. Autoantibody profiles in the sera of European patients with myositis. *Ann Rheum Dis* 2001; 60(2):116-123.
- Castro C, Gourley M. Diagnostic testing and interpretation of tests for autoimmunity. *J Allergy Clin Immunol* 2010; 125(2 Suppl 2):S238-S247.
- Defendenti C, Atzeni F, Spina MF, Grosso S, Cereda A, Guercilena G et al. Clinical and laboratory aspects of Ro/SSA-52 autoantibodies. *Autoimmun Rev* 2011; 10(3):150-154.
- Eriksson C, Kokkonen H, Johansson M, Hallmans G, Wadell G, Rantapaa-Dahlqvist S. Autoantibodies predate the onset of Systemic Lupus Erythematosus in northern Sweden. *Arthritis Research & Therapy* 2011; 13(1):R30.
- Haugbro K, Nossent JC, Winkler T, Figenschau Y, Rekvig OP. Anti-dsDNA antibodies and disease classification in antinuclear antibody positive patients: the role of analytical diversity. *Ann Rheum Dis JID* - 0372355 2004; 63(4):386-394.
- Ippolito A, Wallace DJ, Gladman D, Fortin PR, Urowitz M, Werth V et al. Autoantibodies in systemic lupus erythematosus: comparison of historical and current assessment of seropositivity. *Lupus* 2011; 20(3):250-255.
- Isenberg DA, Manson JJ, Ehrenstein MR, Rahman A. Fifty years of anti-ds DNA antibodies: are we approaching journey's end? *Rheumatology (Oxford)* 2007; 46(7):1052-1056.
- Kattah NH, Kattah MG, Utz PJ. The U1-snRNP complex: structural properties relating to autoimmune pathogenesis in rheumatic diseases. *Immunol Rev* 2010; 233(1):126-145.
- Kumar Y, Bhatia A, Minz RW. Antinuclear antibodies and their detection methods in diagnosis of connective tissue diseases: a journey revisited. *Diagn Pathol* 2009; 4:1.
- Meroni PL, Schur PH. ANA screening: an old test with new recommendations. *Ann Rheum Dis* 2010; 69:1420-1422.
- Petri M, Magder L. Classification criteria for systemic lupus erythematosus: a review. *Lupus* 2004; 13(11):829-837.
- Poole BD, Schneider RI, Guthridge JM, Velte CA, Reichlin M, Harley JB et al. Early targets of nuclear RNP humoral autoimmunity in human systemic lupus erythematosus. *Arthritis Rheum* 2009; 60(3):848-859.
- Putova I, Dostal C, Becvar R. Prevalence of antinucleosome antibodies by enzyme-linked immunosorbent assays in patients with systemic lupus erythematosus and other autoimmune systemic diseases. *Ann N Y Acad Sci* 2007; 1109:275-286.
- Reveille JD. Predictive value of autoantibodies for activity of systemic lupus erythematosus. *Lupus JID* - 9204265 2004; 13(5):290-297.
- Simon JA, Cabiedes J, Ortiz E, Alcocer-Varela J, Sanchez-Guerrero J. Anti-nucleosome antibodies in patients with systemic lupus erythematosus of recent onset. Potential utility as a diagnostic tool and disease activity marker. *Rheumatology (Oxford)* 2004; 43(2):220-224.
- Sinclair D, Saas M, Williams D, Hart M, Goswami R. Can an ELISA replace immunofluorescence for the detection of anti-nuclear antibodies?—The routine use of anti-nuclear antibody screening ELISAs. *Clin Lab* 2007; 53(3-4):183-191.
- Tozzoli R, Bizzaro N, Tonutti E, Villalta D, Bassetti D, Manoni F et al. Guidelines for the laboratory use of autoantibody tests in the diagnosis and monitoring of autoimmune rheumatic diseases. *Am J Clin Pathol* 2002; 117(2):316-324.
- Maidhof W., Hillas O. Lupus: an overview of the disease and management options. *P T* 2012; 37(4):240-9.
- Hahn BH, McMahon MA, Wilkinson A, Wallace WD, Daikh DI, Fitzgerald JD et al. American College of Rheumatology guidelines for screening, treatment, and management of lupus nephritis. *Arthritis Care Res (Hoboken)* 2012; 64(6):797-808.


1

2

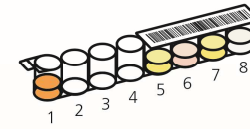3

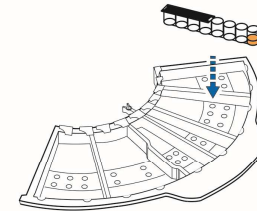4