

Dräger DrugTest® 5000 Drug Testing System for Law Enforcement

The Dräger DrugTest® 5000 system is a fast, simple and accurate way of detecting harmful drugs through oral fluid sampling. Our device has been selected as the only approved drug detection analyzer for law enforcement use in Canada, capable of detecting cocaine and cannabinoids. The analyzer offers easy data management with Dräger Diagnostics Software.



Benefits

Easy Evaluation

The DrugTest 5000 system is designed for easy drug screening and accurate results.

- Quick and simple to operate
 - User-friendly design
 - On-screen instructions guide you through the drug test
-

Quick Confirmation

The device controls analysis time and internal temperature for accurate results. Positive results can be confirmed by a blood sample via laboratory analysis.

Two Step System

System components include:

- DrugTest 5000 Test Kit: a ready-to-use test cassette with a built-in sample collector and volume adequacy indicator.
 - DrugTest 5000 Analyzer: a rechargeable, fully automated analyzer providing controlled accurate analysis on site.
-

Stores, Manages, and Prints Results

The analyzer stores the last 500 test results, along with the date and time. You can enter additional subject or test details using the compact keyboard. Data can be downloaded to a PC or printed using the Mobile Printer.

Non-Invasive

To protect both you and the donor, sample collection is completely non-invasive. When the cassette is handled correctly, you will not come into contact with the donor's oral fluid.

Accessories



ST-3451-2003

Compact keyboards

- with PS2 interface, dimensions approx. 28.2 x 13.2 x 2.4 cm.
- "QWERTY" - English layout, REF: 8315497
- "AZERTY" - French layout, REF: 8315142



ST-8404-2006

Dräger Mobile Printer

The Dräger Mobile Printer prints breath alcohol and drug test measurement results directly on the spot.



D-52008-2012

Carrying case for Dräger DrugTest® 5000

The convenient carrying case for the easy transport of the DrugTest® 5000 system provides space for all required materials: DrugTest® 5000 Analyzer, Mobile Printer, keyboard, test kits and accessories.




D-4337-2014

Hard-shell carrying case

For the secure transport and shipping of the Dräger DrugTest® 5000 e.g. for periodic calibration. Offers storage for printer paper, mouthpieces, keyboard and power cable.

Services

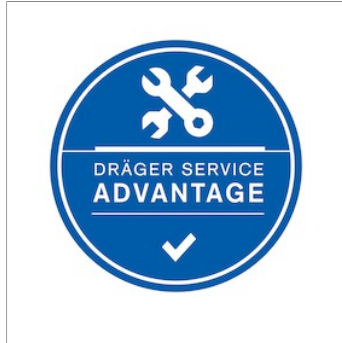

Dräger Service

When you rely on Dräger breath alcohol and drug screening equipment, you can rest assured you made the right choice. For more than 60 years, we've been the global leader in advanced breath alcohol testing and employ cutting-edge technology to provide devices that are fast, reliable, and easy to use. Plus, our devices are backed by world-class service—that's the Dräger Service Advantage.

Related Products



Dräger Alcotest® 6820

Keeping impaired drivers off the road is a tough, unpredictable job. Dräger designed the Alcotest® 6820 specifically to meet the stringent requirements of professional breath-alcohol analysis. The device is fast, easy to use, highly reliable, and approved by the Canadian government as an official screening device for use in law enforcement.



Dräger Alcotest® 5000

The Dräger Alcotest® 5000 is a professional breath-alcohol tester which detects the presence of alcohol. The high-speed breathalyser lets you perform numerous tests in minimal time. Its special funnel reduces the back-flow of expired air to a minimum, preventing the risk of infection to subsequent test subjects.

Related Products



D-43786-2015

Dräger Alcotest® 9510

The Dräger Alcotest® 9510 is an advanced breath alcohol measuring instrument for evidential applications. It is designed to comply with national and international requirements and regulations (OIML R 126). With dual sensor technology, intuitive colour touch screen interface, and intuitive design, the Alcotest® 9510 sets superior standards.



D-6296-2009

Dräger SSK 5000

Ideal for surface drug testing applications, the Dräger SSK 5000 sampling system is quick and easy to use. When used with the Dräger DrugTest 5000, the sample can be transferred to a test cassette and screened for possible drug contamination. The sample can also be stored and transported for further testing.

Technical Data

Supply Voltage	12 V DC (11 to 15 V DC)
Current Consumption	Typically 3 A
Dimensions (W x H x D)	7.9" x 9.8" x 8.7"
Weight	9.9 lbs
Temperature	Operation: 5 to 40 °C (+41 to +104°F) Storage/Transportation: -20 to +60°C (-4 to +140°F)
Moisture Range	5 to 95%rh, non condensing
Interfaces	Optical IRDA (printer), PS/2, USB Slave
Duration of One Measurement	Between 4 to 5 minutes
Storage Capacity	500 data sets
Carrying Bag	For Mobile System

Ordering Information

Item	Order No.
Dräger DrugTest® 5000 Analyzer including power supply, connecting cable and instructions for use.	37 01 413
Consumables	
Dräger DrugTest® 5000 Test Kits - 20/Box (COC, THC)	37 00 791
Dräger DrugTest® Surface Screening Kits - 2 Panel -20/Box	83 20 490
Accessories	
Compact Keypad (PS/2, dimensions approx. 11.1" x 5.2" x 0.9") "QWERTY" English Keyboard Layout	83 15 497
Car Cable 12 V (Connecting Cable from car socket/cigarette lighter to DrugTest® 5000 Analyzer)	83 12 166
Dräger Mobile Printer	83 19 310
Thermal Paper for Dräger Mobile Printer (5 roll)	83 19 002
Carrying Bag for Dräger DrugTest® 5000 Analyzer	83 22 675
USB Connection Cable for Communication with a PC	AG02661
Dräger DrugTest® Training Test Kit (TTK) - 6 pcs.	83 23 640

An optional carrying bag is available for the safe transportation of all system components. The keyboard enables information on the test subject and system operator to be entered, allowing a record of this data to be printed out on the spot along with the results.

Notes

Notes

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Allemagne
www.draeger.com

CANADA

Draeger Safety Canada Ltd.
2425 Skymark Avenue, Unit 1
Mississauga, Ontario, L4W 4Y6
Tél +1 905 212 6600
Toll-free +1 877 Drager 1
(+1 877 372 4371)
Fax +1 877 651 0902
Toll-free Fax +1 800 329 8823

Locate your Regional
Sales Representative at:
www.draeger.com/contact

